



Photo: is courtesy of Enable-Me

COMFORT • SUPERIOR SEATING • ACCESS • ADJUSTABILITY

REHAB SHOWER COMMUNE CHAIRS
The Choice of Health Care Professionals

RAZ-CAT COMPACT ATTENDANT TILT



- > Visco Foam Seat, 18"W with Oblong Aperture
- > Adjustable Seat Depth: 11" - 16"
- > Back Widths Available: 14", 16", 18"



- > Z7MFX12N Adjustable 12"-16.5" Footrests (shown) Z7MFX8N Adjustable 8"-12.5" Footrests (not shown) Z708ADJ Adjustable Heel Loops (optional)



- > Small Lateral Support with Depth Adjustable Pads Offset is Adjustable between 1"-3"
- > Item Code: Z740AP1 L/R



- > Large Lateral Support with Depth Adjustable Pads Offset is Adjustable between 1"-3"
- > Item Code: Z7 40AP2 L/R



- > Small Hip Adductor with Depth Adjustable Pads Offset is Adjustable between 0"-2.5"
- > Item Code: Z747AP1 L/R



- > Large Hip Adductor with Depth Adjustable Pads Offset is Adjustable between 0"-2.5"
- > Item Code: Z747AP2 L/R



RAZ has created the most configurable, adjustable and comfortable rehab shower chairs on the market today. The chairs were developed with meticulous attention to detail — from the electro-polished stainless steel frames with easy access for users and caregivers, to precisely contoured seats molded with soft integral skin foam. Only RAZ offers a commode seat that adjusts to accommodate pelvic rotation with a centering of the aperture around the user's ischial tuberosities.

RAZ Rehab Shower Commode Chairs, provide positioning, pressure reduction, promote function and are ideal for the rehab client and the institutional user.



RAZ-AT

WHEN SEATING IS IMPORTANT, CHOOSE RAZ...

ISCHIAL PELVIC ALIGNMENT SEAT (IPAS)

Recognizing that many users spend lengthy periods of time sitting, we knew the seating for the RAZ Rehab Shower Commode series was a critical element. The IPAS was developed with contours to increase pelvic stability and reduce pressures over the ischials and coccyx. The IPAS provides 2" of seating depth adjustment and can be fitted with optional seat spacers for seat angle adjustment. The IPAS is available in 16" or 19" depths with open or bridged front.


IPAS BRIDGED FRONT SEAT

IPAS OPEN FRONT SEAT



PRESSURE REDUCTION SEAT (PRx)

For users with more complex seating requirements. The PRx is made with a high density foam base foundation, designed to provide superior stability and contours, while the surface layer of self-skinning gel foam provides a soft interface against the body. The PRx seat is available in ten configurations including left, right and rear cut-out for improved access or full seat models for comfort and pressure reduction, all in 16" or 19" depths.

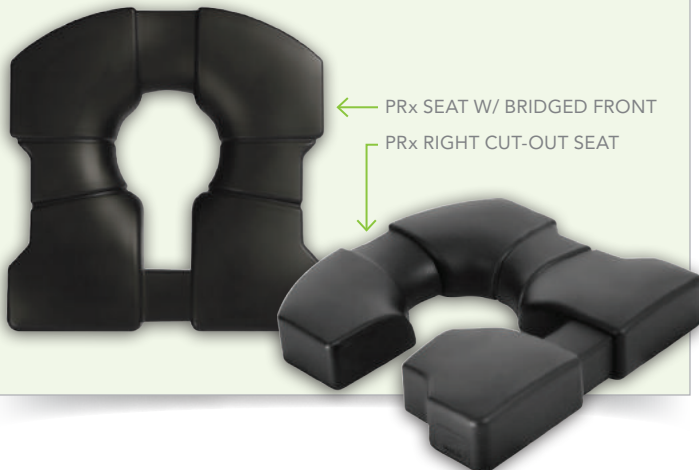

PRx SEAT W/ BRIDGED FRONT

PRx RIGHT CUT-OUT SEAT

Exclusive RAZ-AT Features


ADJUSTABLE PIVOT POSTS

- > Provides 4" of seat height adjustment without frame-toilet bowl interface
- > Pressurized gas springs work against gravity making tilting up to 40° easier for an attendant



TILT ASSIST PEDAL

- > Enables attendant to use additional strength from legs to tilt larger/heavier users

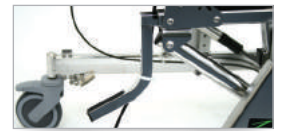


NEW! Z515 RECLINING BACK

- > 30° of recline
- > 20" back height
- > 25° of tilt




SEAT SPECIFICATIONS

01. Seat width
02. Seat depth
03. Aperture width
Standard seat: 7.25"
PRx seat: 6.25"
04. Aperture length
Standard seat: 7"
PRx seat: 6"
05. Front opening
Standard seat: 3.5"
PRx seat: 2.75"
06. Optional front access
07. Optional 3" bridge
08. 1" side cut-outs

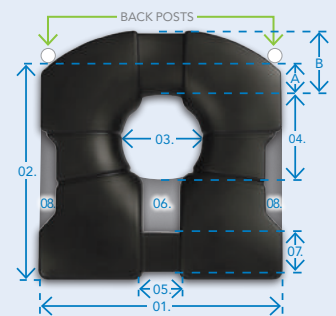

PRx side & rear opening

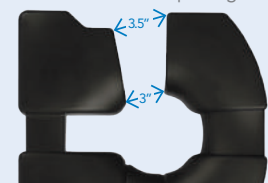

RAZ AT, SP & AP
COME STANDARD WITH THE
HARMONY
BACKREST
VINYL UPHOLSTERY IS
A NO-CHARGE OPTION



RAZ-SP



RAZ-AP



NEW! RAZ-AT600

Exclusive RAZ-AP/SP Features



ADJUSTA-BACK

- > Adjusts to 85°, 90°, & 95° back angle to suite user's needs
- > Armrest adjusts in height/angle
- > Includes locking armrest feature
- > Available in 17" & 20" back height



MULTI-POSITION AXLE PLATES

- > Provides up to 4" of vertical seat height and 3.5" of horizontal adjustment for optimal user-toilet bowl access

RAZ-AT600 Features

- > Weight assist spring system
- > 25° of tilt
- > Seat widths from 22" - 30" wide
- > Optional bridged or open front or side and rear opening

RAZ REHAB SHOWER COMMUNE CHAIR FEATURES



01



02



05



06




07



08



09

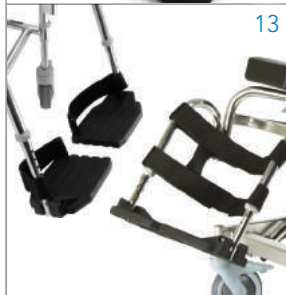

13




14



15




16

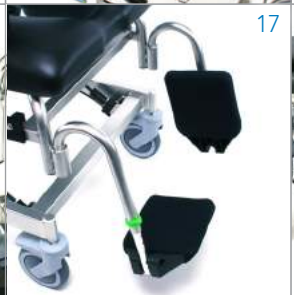

17




Exclusive RAZ-AP600/SP600/AT600 Features

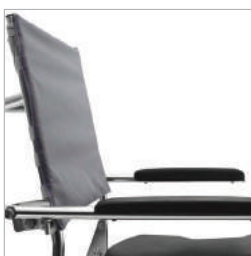

- MULTI-POSITION BACK**
- > Adjusts to 90° & 100° back angle for personalized comfort
 - > Flip up/fixed height armrests
 - > Armrests adjust in angle based on user's unique needs
 - > Available with 20" back height



- CENTRE OF GRAVITY ADJUSTABLE FRAME**
- > Provides independent horizontal adjustment of back (1.5") and axle plate (1.5") for optimal maneuverability and stability
 - > Axle plates provide up to 3.75" of seat height adjustment

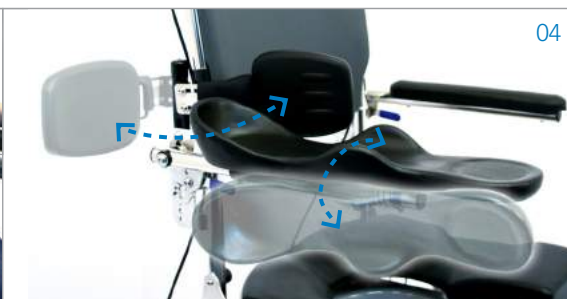

- 01. Harmony Back
- 02. Adjustable Tension Vinyl Backrest
 - N/C option. Spray or wipe for easy cleaning
- 03. Body Point Lap & Chest Belt
- 04. Swing Away Lateral Support & Pivoting Arm Mount
 - Armrest pivots down to allow armrest to flip up
- 05. Flip Back 1-pc Foot Rest
- 06. V Style 1-pc Foot Rest
- 07. E & J Replica Seat
- 08. Visco Interface Seat
 - with optional bridged or open front or right, left or rear opening
- 09. Leg & Hip Adductors
 - Height, depth & length adjustable
- 10. Adjustable Seat
 - Provides 2" of depth adjustment
- 11. Commode Pail
 - Fits most Autoclaves
- 12. Molded Flat Arm Pad
 - 3 mounting tracks allow multiple mounting positions
- 13. Heel Loops/H Straps
- 14. Body Point Calf Panel & Calf Strap
- 15. Depth & Angle Adjustable Foot Plate
- 16. MFX 8 & 12 Foot Rests
 - MFX 8: 85° angle & adjusts from 8"-12" in length
 - MFX 12: adjusts 80° & 87° & from 12"-16" in length
- 17. Neoprene Footplate Covers
- 18. Elevating Leg Rests (FOR 1819 & 1816 SEATS)
- 19. AP & AT Directorial Caster Locks
- 20. Large Head Rest Pad 12" x 6"

SPECIFICATIONS

STANDARD CONFIGURATION	RAZ-AP	RAZ-SP 22" 24" wheels	RAZ-AT	RAZ-AT600	RAZ-AP600	RAZ-SP600	RAZ-CAT
BACK HEIGHT <i>(from commode seat)</i>	17" standard 17" optional 20" optional		17" standard 20" optional 24" optional	21"	20"		15" & 18"
SEAT HEIGHT <i>(front seat to floor)</i>	19.5" - 23.5"	19.5" - 22.5" 19.5" - 23.5"	21" - 25"	22.5" - 23.5"	5": 19.25" - 23.25" 6": 20.25" - 24.25"		21" - 25"
SEAT WIDTH	18"			23"		18"	
SEAT DEPTH	16" or 19"			18"		11" 12" 13" 14" 15" 16"	
SEAT ADJUSTMENT <i>(fore and aft)</i>	Up to 2"			Up to 1.5"		-	
OVERALL WIDTH	23" w/Std. Arm Pad 25.5" w/Z715 Arm Pad	26" <i>(outside hand rim)</i>	23" w/Std. Arm Pad 25.5" w/Z715 Arm Pad	24.75"	27" w/Std. Arm Pad 29.5" w/Z715 Arm Pad	30" <i>(outside hand rim)</i>	21.8" 23"
CASTER SIZE & TYPE	Four 5" dual locking	Two 5" dual locking	Four 5" dual locking	Four 6" dual-locking	Four 5" dual locking Four 6" dual locking	Two 5" dual locking Two 6" dual locking	Dual-locking 5" casters (x4)
REAR WHEELS <i>(treaded urethane tires)</i>	-	22" standard 24" optional	-	Optional retrofit kit		24"	-
OVERALL LENGTH <i>with foot rests</i>	40.2"	40.3" - 43.8" 41.1" - 44.6"	45.25"	42" - 43"	39" - 40.5"	44" - 45.25"	44.4"
OVERALL LENGTH <i>without foot rests</i>	28.6"	28.7" - 32.2" 29.5" - 33.0"	33.6"	32.5"	29.5"	32" - 34.5"	31.9"
FOOT REST ADJUSTMENT <i>(from seat to footplate)</i>	Z705N: 16" - 20.5" ZFREXT: 18" - 22.5"			14.5" - 19" 14.5" - 19"		12" - 16" 14" - 18"	
WIDTH BETWEEN ARMRESTS	18.5"			22.75"	22.5" 22.5"		14.75" 16.75" 18.75"
ARMREST HEIGHT	Z505 back: 8.5" Z506 back: 9.5" Z510 back: 7" - 9.5" Z515 back: 8" - 10.5"		Z505 back: 8.5" Z506 back: 9.5" Z510 back: 7" - 9.5" Z515 back: 8" - 10.5" Z520 back: 8.5" - 10.5"	Adjustable: 9" - 10"	8"		6" - 8" 7" - 9"
CLEARANCE HEIGHT UNDER SEAT	14.8" - 18.8" w/pan bracket 16" - 20" w/o pan bracket	10" - 13" w/pan 14" - 17" w/o pan 15" - 18" w/o commode pan bracket	14.8" - 18.8" w/pan bracket & 22" wheel 16" - 20" w/o pan bracket & 22" wheel 14.8" - 17.8" w/pan bracket & 24" wheel 16" - 19" w/o pan bracket & 24" wheel	18" - 19" w/pan 18.25" - 19.25" w/o	5": 14" - 18" w/pan 5": 15.5" - 19.5" w/o 6": 15" - 19" w/pan 6": 16.5" - 20.5" w/o	5": 15" - 18" w/pan 5": 16.25" - 19" w/o 6": 15" - 18" w/pan 6": 16.5" - 19.5" w/o	-
WEIGHT CAPACITY	300 lbs			600 lbs		250 lbs	
PRODUCT WEIGHT	38 lbs	47 lbs	65 lbs	85 lbs	5": 54 lbs 6": 56 lbs	5": 68 lbs 6": 69 lbs	64 lbs.
FRAME MATERIAL & HARDWARE	Stainless Steel						

WARRANTY Two year limited warranty on the frame against manufacturer's defects.
One year on remaining components subject to normal wear and tear exclusions.



Dynamic Health Care Solutions
753011 Second Line EHS
Orangeville, ON L9W 5W4
Tel: (866) 875-2877 | (519) 942-8441
Fax: (866) 875-2878 | (519) 942-2186
DynamicHCS.com